




Features of Magnetic Infrared Warm Heater NM-70

- It relieves pain in knee and shoulder by simple handling.
- Temperature accommodation is possible from 40°C to 70°C.
- It has human engineering design to protect knee.
- It helps blood circulation with built-in magnetic of 800 gauss for medical use.
- It maximizes infrared effect with tourmaline ceramic.





JOINT CARE NM-70

Manual for NM-70(JOINT CARE)

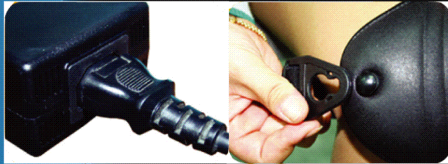


1. Connect power cable and wear the device on the section to be used, using fixation band.

2. Connect power jack in adaptor to remote control cable.

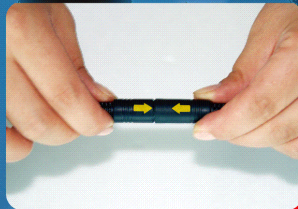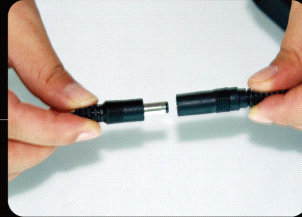

3. Connect warmer cable to remote control cable.
Note: Connect them to arrow direction of each other.

4. Set your desired stimulation time after pushing power button in remote controller.
(Setting time is available for 60 min., 90 min., 120 min., and 180 min. It is automatically off when reaching setting time.)

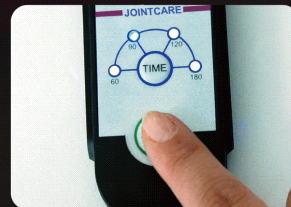

5. Set your desired stimulation temperature.
(Setting temperature is available for 40~65°C. The unit is increasing by 5°C under 50°C and increasing by 2°C over 50°C.)
For user's feasibility, current temperature is indicated after displaying setting temperature at temperature display window for 3 seconds when increasing temperature.
However, in lowering temperature, setting temperature is displayed only.



6. Setting of time and temperature is available even in the middle of operating.
Pushing POWER in operating, the movement stops.

Note: Use the device after wearing it on the protection supplied by us or outer garment.



Caution to use

1. If any trouble occurs, contact skilled technician after indicating and recording the state of trouble, not handling the heater, and record.
2. Do not alter or remodel the heater.
3. Avoid to use the heater in the place where big current is used at a moment such as generator and X-ray equipment or electronic waves are largely exposed.
4. Handle it by fixed order and put Power switch Off after using.
5. Do not use out of purpose.
6. Keep regulated stimulation time.
7. Do not apply it to the person who uses artificial heart device (pace maker).



Specification

1. Medical product name: Warm heater for individual
2. Manufacturer: 1642-10, Donghwa Newest Medical Machinery Industrial Site, Donghwa 3ri, Munmak-eup, Wonju-si, Gangwondo
3. Certificate of manufacturer: 1478
4. Item permission (Application) No.: 06-873
5. Weight: 2.6Kg, 1set
6. Rating voltage and frequency:
DC 24V / 2.5A (DC Adaptor Input voltage and rating frequency: AC 100-240, 50-60Hz)
7. Electricity consumption in Max. output: Under 60W
8. Protection Type to electric shock: Grade II machine, B type.