



CENTRE FOR TESTING AND CERTIFICATION - MECH-TEST

Mechanical Laboratory

05-077 Warszawa-Wesoła, ul. Klonowa 22

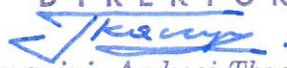
tel.: +48 603 23-26-45, e-mail: cbc.mech.test@gmail.com, www.cbc.org.pl

Date 16.06.2021

TEST REPORT NO. **CBC-145/2021**

Page 1 of 4

Subject of testing:	<i>Walking aids with built-in handgrips and three or more legs of which two or more are having wheels, which provide support whilst walking</i>	Classification according to PN-EN ISO 9999:2017-02 : 12 06 06
Type / Model:	<i>Leopard rollator with soft wheels, anthracite grey, seat height 62cm, incl. backbelt and cane holder</i> <i>SOFT rear wheel with black colored rim for rollator (1 pcs.) with RolloGuard brake system</i>	Item no.: 312306 Item no.: 31233R SN: 0205
Manufacturer:	<i>MOBILEX A/S Grønlandsvej 5 DK-8660 Skanderborg</i>	Number of specimens: 1
Applicant:	<i>A-Net s.c. 93-469 Łódź, ul. Łaskowice174</i>	
Kind of testing	<i>Mechanical testing for conformity with PN-EN ISO 11199-2 : 2005</i>	
Test started:	<i>8.06.2021</i>	
Test finished:	<i>16.06.2021</i>	

Approved by:
DYREKTOR

mgr inż. Andrzej Tkaczyk

Special comments / enclosures:

Copyright © 2012 by Centre for Testing and Certification (applicable to report form)

Test results refer only to tested units.

Test results reported here are not applicable to the further modifications of the product affecting its structure, material or technology.

This test report shall be neither copied differently as in the whole nor be published without written consent of the Laboratory.



CHARACTERISTIC OF PRODUCT

Name : <i>Leopard rollator</i>	Dimension of rollator: --
Item no: <i>312306</i>	SN: <i>0205</i>
Maximum permissible user mass: <i>150 kg</i>	Mass of rollator: <i>9,04 kg</i>

Description		Comments	
Elements/parameters/materials/dimensions			
Dimensions of walking rollator (fig. 2 PN-EN ISO 11199-2)	Distance between handgrips (dimension 2)	<i>475-470 mm</i>	
	Angle between of handgrip axis and direction of movement (α)	<i>0°</i>	
	Height of rollator (dimension 6)	<i>767/740 mm</i>	<i>min.</i>
		<i>1010mm</i>	<i>max.</i>
	Width of rollator (dimension 5)	<i>608 mm</i>	
	Turning width (dimension 1)	<i>845 mm</i>	
	Length of rollator (dimension 4)	<i>664/684 mm</i>	
Dimensions of folded rollator (mm)		<i>664 x 790 x 236</i>	
Fig. 3	Handgrip - diameter	<i>38 mm</i>	<i>Anatomical handgrip</i>
	Handgrip - length	<i>120 mm</i>	
Wheels of rollator	Front wheels – quantity	<i>2</i>	<i>castor wheels</i>
	Front wheels - diameter	<i>200 mm</i>	
	Front wheels – width	<i>31/36 mm</i>	
	Front wheels - brake	<i>none</i>	
	Rear wheels – quantity	<i>2</i>	
	Rear wheels - diameter	<i>200 mm</i>	
	Rear wheels - width	<i>31/36 mm</i>	
	Rear wheels - brake	<i>Included</i>	
Tip	Diameter		
	Material	<i>Not any</i>	
	Colour		
Material of rollator (fig. 1)	Front legs	<i>Aluminum,</i>	
	Bracing member (no. 8)	<i>Hard plastic,</i>	
	Rear legs	<i>Bolts, nuts</i>	
	Height adjusting device (no. 4)		
	Handgrip (no 5), Brake elements	<i>Hard plastic</i>	



Mobilex Leopard rollator
 a better life in motion For outdoor and indoor use
 Item no. 312306

MOBILEX A/S Max. 150kg
 Grønlandsvej 5
 DK - 8660 Skanderborg Max. width = 60.0cm
 Tel: +45 87 93 22 20 Max. load in basket = 5 kg
 www.mobilex.dk

CE Produced 2020.10 Serial no. 0205

SN 3123062020100205



RESULT OF TESTS ACCORDING TO PN-EN ISO 11199-2:2005

Requirements according to clause	Test method according to clause	Checked characteristics/assemblies/parameters		Real value	Test result	Comments	
4.1	Measur.	Manoeuvrability		\varnothing 200 mm width 31mm Conf.	Pos.	\varnothing front wheels ≥ 75 mm outdoor intended rollator: \varnothing front wheels ≥ 180 mm width of wheels ≥ 22 mm	
4.2	5.3	Forward-direction stability		18,9 ⁰ Conf.	Pos.	Stability required $\geq 15^\circ$	
4.2	5.4	Backward-direction stability		9,8 ⁰ Conf.	Pos.	Stability required $\geq 7^\circ$	
4.2	5.5	Sideway-direction stability		4,2 ⁰ Conf.	Pos.	Stability required $\geq 3,5^\circ$	
4.2	5.6	Stability – with loaded basket, bag, drip, oxygen cylinder		forwards	17,2 ⁰ Conf.	Pos.	Stability required $\geq 15^\circ$
				backwards	13,6 ⁰ Conf.	Pos.	Stability required $\geq 7^\circ$
				side	5,4 ⁰ Conf.	Pos.	Stability required $\geq 3,5^\circ$
4.3	V/I	Brakes	Servicing facility during rollator motion with more than 2 wheels		Conf.	Pos.	
	V/I		Parking brakes in rollator with more than 2 wheels and resting seat or intended for outdoor use		Conf.	Pos.	
	5.7.1.1		Brake grip distance (fig. 4, dimension 1)		65 mm Conf.	Pos.	≤ 75 mm
	5.7.1		Running brake effectiveness		Conf.	Pos.	Movement of rollator ≤ 10 mm in 1 minute
	Measur.		Force to set parking brake		30N Conf.	Pos.	≤ 60 N
	Measur.		Force to release parking brake		20N Conf.	Pos.	≤ 40 N
	5.7.2		Parking brake effectiveness		Conf.	Pos.	Movement of rollator ≤ 10 mm in 1 minute
	V/I		Possibility to compensate brake wear		Conf.	Pos.	
V/I	Brake not adversely affected by folding, unfolding or adjusting actions of rollator		Conf.	Pos.			
4.4	Measur. V/I	Handgrip		38 mm Conf.	Pos.	Width of handgrip ≥ 20 mm and ≤ 50 mm	
4.5	Measur. V/I	Leg section and tip		--	N/A	\varnothing tip ≥ 35 mm (tested rollator is equipped in four wheels)	
4.6	5.10	Resting seat – static loading durability		Conf.	Pos.	1 minute under load 1,2 x user's weight $\pm 2\%$ (180kg)	
4.7	5.12	Mechanical durability	Fatigue test		Conf.	Pos.	200 000 cycles with load. 120kg $\pm 2\%$, f=1Hz
4.7	5.11		Static loading test		Conf.	Pos.	loading 180kg $\pm 2\%$, 5sek. NOTE 1
4.8	V/I	Adjusting devices		Conf.	Pos.		
4.9	5.14	Folding mechanism		Conf.	Pos.		
4.11	ISO 10993-1	Materials and finish	Biocompatibility of material with human body		--	N/T	
	V/I		Free of discolouring of skin or clothing in contact with rollator materials		Conf.	Pos.	
	V/I		Burr, shar edges, projections		Conf.	Pos.	
Marking and labelling of product							
6.2	V/I	a) Maximum user mass		--	N/T		
		b) Maximum safe working load (SWL) to be marked on accessories		--	N/T		
		c) Maximum allowed angle between the longitudinal centreline of the handle and the direction of motion, if the handles are sideways adjustable		--	N/T		
		d) Manufacturer's name or trade name and address		--	N/T		
		e) Manufacturer's model identification name and/or number		--	N/T		
		f) Month and year of manufacture		--	N/T		
		g) Maximum extension of the height adjustment, marked on the adjusting members		--	N/T		
		h) Maximum width of the rollator		--	N/T		
		i) Rollator intended for outdoor/indoor use		--	N/T		
4.10	V/I	Warning showing allowed angle between handle axis and direction of movement or physical stop of angle adjusting		--	N/T		

Contents of user manual and/or assembly manual or clear and indelible marking of product

6.3	V/I	a) Maximum rollator height	--	N/T
		b) Minimum rollator height	--	N/T
		c) maintenance and cleaning instructions, including a description of the method and suitable cleaning agents and any precautions needed to avoid corrosion and/or ageing of the materials used in construction of the rollator	--	N/T
		d) Instructions for assembly, adjustment of all kinds, folding and unfolding	--	N/T
		e) Warnings and advice about precautions relating to safe distances between moving and stationary parts (see EN 12182, Clauses 12 and 13, for guidance)	--	N/T
		f) Maximum safe working load (SWL) for load carrying accessories such as basket, tray, shopping bag, etc.	--	N/T
4.10	V/I	Warning in user manual on consequences of such an adjustment of angle between handle longitudinal axis and direction of movement outside allowed value (when handles are adjustable aside).	--	N/T

TEST CONDITIONS

Ambient temperature	19°C	Required temperature 21°C ±5°C
Relative humidity of air:	55 %	Not required
<i>Comments:</i>		
All tests performed with maximum height adjustment of rollator.		
All tests performed in the least stable position of self-adjusting wheels.		
Tests performed with handles positioned at their maximum (allowed) angle to the direction of motion (when adjustment is possible).		
Sequence of tests: stability test, static loading test, fatigue test.		
One rollator was tested.		
During visual inspection before testing any visible defects that could have influence on test results were not stated.		

Pos. – positive; Neg – negative; N/T – not tested; N/A – not applicable; N/R – not required, N/O – not occurred, V/I.- visual inspection, Conf.- conformed.

NOTE 1: Deformation – 30 mm, elastic deformation – 30 mm, permanent deformation – 0 mm (0,0%).

NOTE 2: Also applies to the centrifugal brake.

CONCLUSIONS:

Testing object **conforms** with requirements of PN-EN ISO 11199-2 : 2005, in scope of mechanical testing ordered by client, excluding biocompatibility tests of material with human body according to PN-EN ISO 10993-1:2010.

---- END ----